


Rubber Fab Technologies Group's

# BIOLOGICAL INDICATOR GASKET

The "New" Smart Gasket

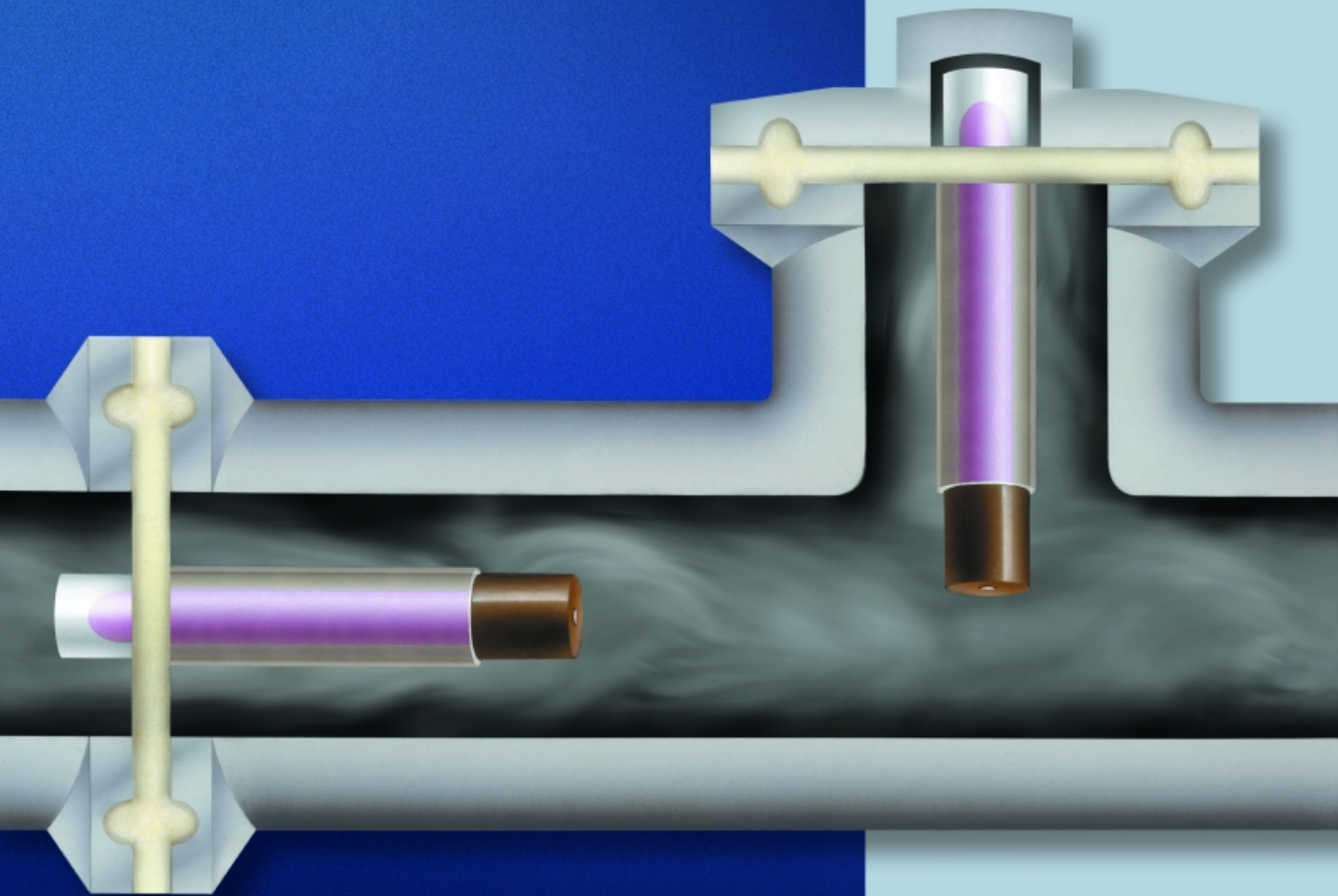

## The Biological Indicator Gasket

Rubber Fab is proud to introduce a new product that enables the use of a **self contained biological indicator** in sanitary process lines. This unique Rubber Fab Smart Gasket® permits either top loading or inline positioning of an **EZTest® self contained biological indicator and/or a temperature probe**.

The great advantage of a **self contained biological indicator** over the traditional spore strip in glassine is **the time saved** when confirming spore kill. **Results are easily obtained in only 48 hours** (versus the five to seven days required with a spore strip) **and without problematic laboratory transfer**.

You can feel secure knowing that the **biological indicator** will not be lost down stream in the process. The fail-safe Rubber Fab Smart Gasket® holds the indicator in place for **easy retrieval**. Available for 3/4", 1.0", 1.5" and 2.0" Sanitary, Tri-Clamp process lines.

## Meets High Pharmaceutical Standards

### Biological Indicator Gasket

- U.S. Pharmacopeia Class VI Certification
- Cytotoxicity Criteria
- CFR Title 21 Section 177.2600
- Traceability: Lot and Batch
- Certification: Lot and Batch
- ASME-BPE Standards
- USDA Standards
- 3-A Sanitary Standards
- Current Good Manufacturing Practices (CGMP)

### EZTest®

- U.S. Pharmacopeia
- ISO 11138
- European Norm 866



Rubber Fab Technologies Group  
Limecrest Road

P.O. Box 626 • Andover, NJ 07821  
phone: 973-579-2959 fax: 973-579-7275

[www.rubberfab.com](http://www.rubberfab.com)



**Additional technical and usage information about the EZTest® biological indicator and related testing products is available at [www.sgmbiotech.com](http://www.sgmbiotech.com).**

## Validation of Steam Systems

- Fermenters
- Pharmaceutical Tanks
- Sanitary Inline Pipe
- Pharmaceutical Hose



## Smart Gasket Technology



The Biological Indicator Gasket is also available with one (1) Smart Gasket Thermocouple Port.

See Rubber Fab's Smart Gasket literature on our website for more information on Rubber Fab's other inline system validation products.

Rubber Fab's Biological Indicator Gasket is patent pending.

Smart Gasket® is a registered trademark of Rubber Fab Technologies Group.

EZTest® is a registered trademark of SGM Biotech, Inc.

**SCM BIOTECH, INC.**

Rubber Fab Technologies Group is a member of:

Distributed by:

